


Hill-Rom NP100 Bidensity

Efficient and safe pressure ulcer prevention

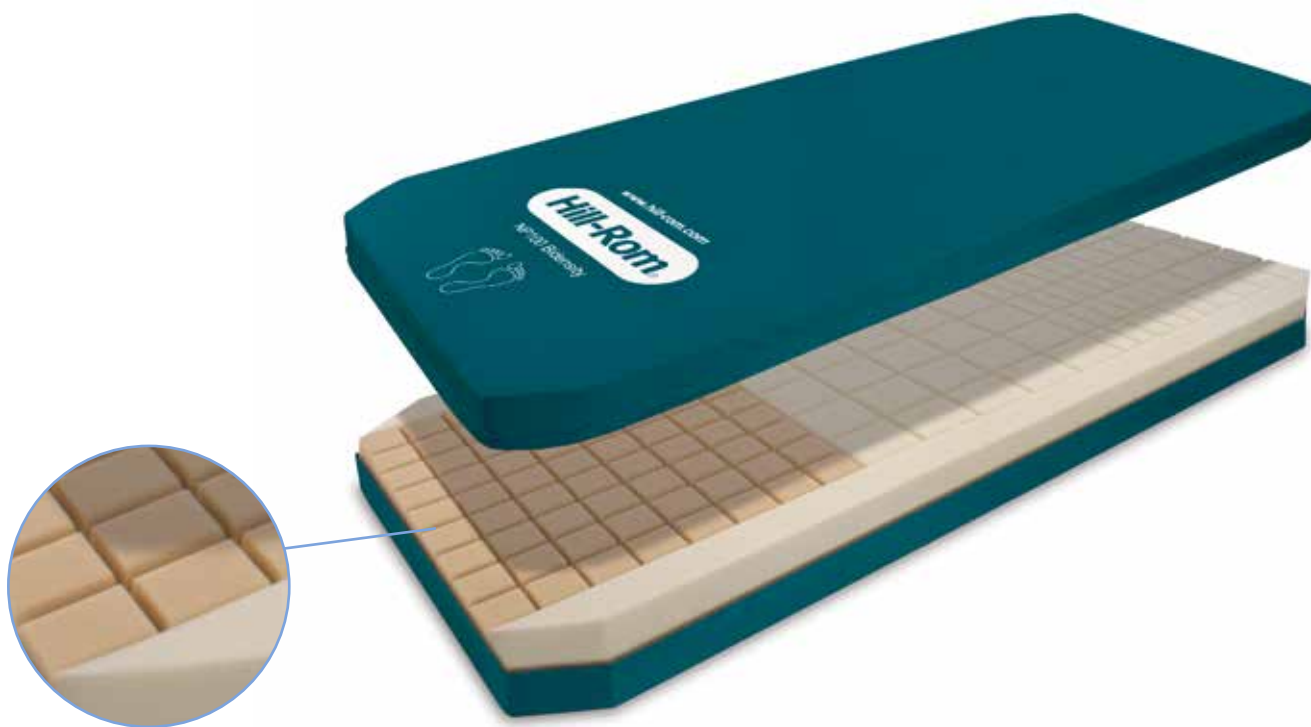

Features and benefits

- Cube-cut design conforms independently to the body shape and reduces interface pressure and friction
- The independent cubes promote airflow to reduce moisture build-up which in turn reduces maceration and shear
- Low-density heel zone offers gentle support for heels and ankles
- Easy to clean and hygienic
- High density side sections designed for patient safety and improved egress
- Adapts to most bed frames and conforms to all bed positions
- Welded cover with built-in carry handles
- X-ray translucent and latex-free

Enhancing outcomes for patients and their caregivers:

Hill-Rom

Hill-Rom NP100 Bidensity

High-resilience, bidensity foam mattress

With features such as independent, bidensity, cube-cut foam, the Hill-Rom NP100 Bidensity mattress reduces interface pressure and friction by individually conforming to the patient's shape. The cube-cut foam promotes airflow under the surface to reduce moisture build up which in turn reduces maceration and shear.

The Hill-Rom NP100 Bidensity mattress is a cost-effective solution to help prevent skin breakdown and pressure ulcers. With its low-density heel zone it also offers gentle support for heels and ankles. It is adaptable to most bed frames and conforms to all bed positions.

The fully enclosed mattress, with its waterproof, yet vapour permeable, multi-stretch, bacteriostatic and fungostatic polyurethane cover, is easy to clean and disinfect enhancing hygiene for patients and caregivers.



The Hill-Rom NP100 Bidensity offers 1 extended length

Standard size: 198 cm

+5 cm

Technical Specifications

Mattress (customised sizes available)

Length	198 or 203 cm
Width	85 or 90 cm
Height	14 cm

Warranty

Hill-Rom offers a 5 year warranty on the foam and a 2 year warranty on the cover. For details of the warranty provided please contact your local Hill-Rom representative.

Safe Working Load

(Maximum Patient Weight)	250 kg
--------------------------	--------

Seams	Welded
-------	--------

Disinfection

Chlorine-based solutions can be used. Concentration must be less than or equal to 1,000 ppm.

Hill-Rom is a leading global medical technology company with more than 7,000 employees in over 100 countries. We partner with health care providers by focusing on patient care solutions that improve clinical and economic outcomes in five core areas: **Advancing Mobility, Wound Care and Prevention, Clinical Workflow, Surgical Safety and Efficiency, and Respiratory Health.** Hill-Rom people, programs, and product brands work towards one mission: **Every day, around the world, we enhance outcomes for our patients and their caregivers.**

Not all products/options are available in all countries. For further information about this product or a service, please contact your local Hill-Rom representative or visit our webpage:

www.hill-rom.com

This document is destined solely for use by healthcare professionals. Medical devices shown in this brochure are intended for use with patients in departments of healthcare establishments.



These medical devices bear the CE mark and comply with all applicable regulations. Hill-Rom recommends you carefully read the detailed instructions for safe and proper use included in the documents accompanying the medical devices. The personnel of the healthcare establishments are responsible for the proper use and maintenance of these medical devices.

Medical beds/mattresses manufacturer:

Hill-Rom SAS, B.P. 14 - Z.I du Talhouët, 56330 Pluvigner - France (also European Representative)
Hill-Rom Industries SA, 188 Rue du Caducée, Parc EuroMédecine, 34195 Montpellier Cedex 5 - France
Völker GmbH, Wullener Feld 79, 58454 Witten - Germany
Hill-Rom Holdings, Inc., 1069 State Route 46 East, Batesville, IN 47006 - USA
Lifts manufacturer: Liko AB, Nedre vägen 100, 975 92 Luleå - Sweden

Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products.

©2015 Hill-Rom Services, Inc. ALL RIGHTS RESERVED.

Doc. No: 5EN116304-02, 28 January 2015

Enhancing outcomes for patients and their caregivers:

Hill-Rom